


Quality built, affordably priced



*Convenient handles for easy movement and positioning.*

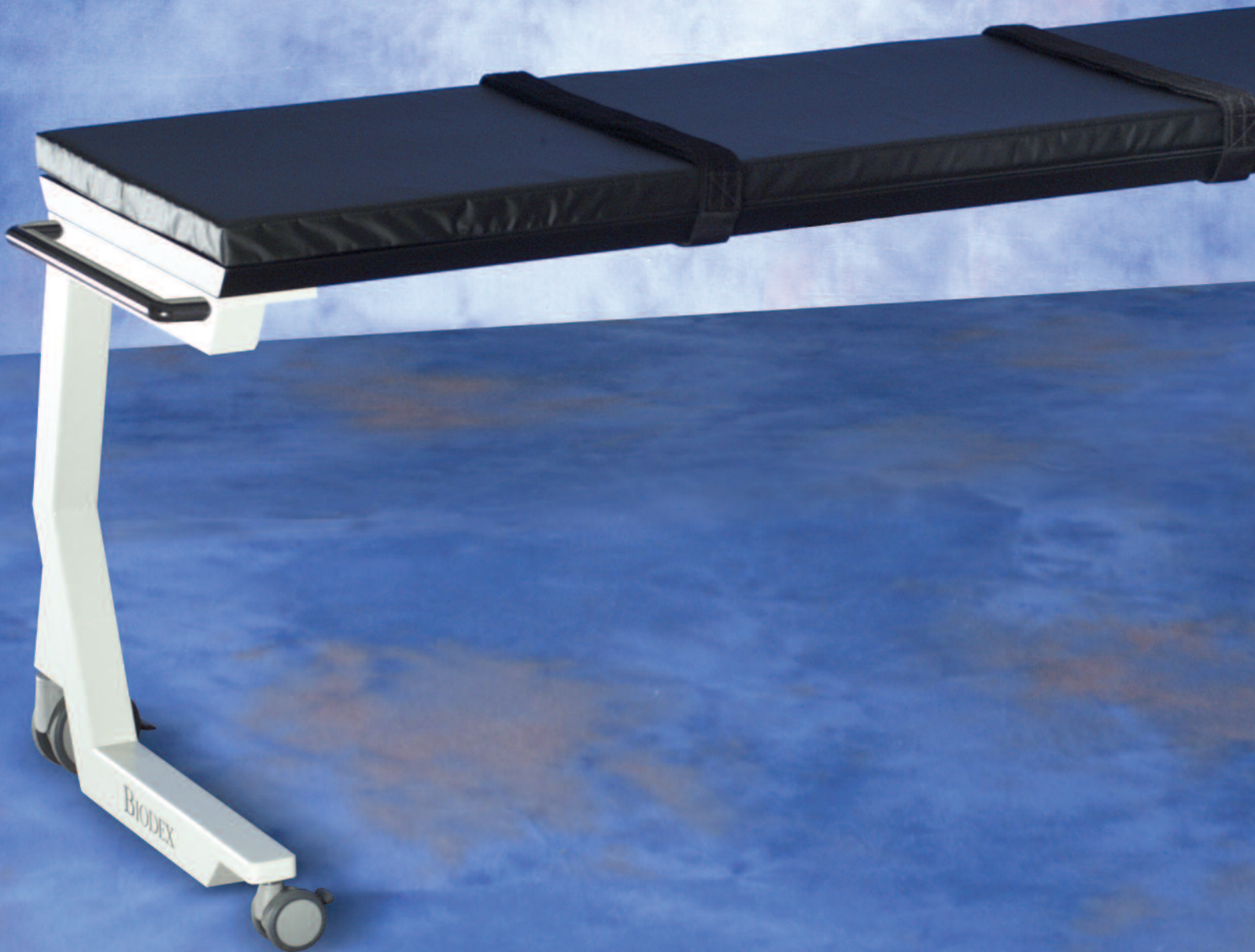

[www.biomed.com/fixedheighttable](http://www.biomed.com/fixedheighttable)

# Unrestricted head-to-toe C-Arm positioning during image-guided procedures

- Extra-large low-attenuation carbon fiber tabletop
- Unrestricted end-to-end C-Arm positioning
- Four twin-wheel locking casters

Designed for unrestricted head-to-toe C-Arm positioning during image-guided procedures, this unit offers structural stability and manufacturing quality at an affordable price. The extra-large 80" x 24" carbon fiber tabletop has no cross members between the two C-shaped supports at each end of the table. This open side-cantilever design provides unobstructed end-to-end C-Arm placement for full head-to-toe fluoroscopic visualization of any size patient.

Two-inch thick table padding assures patient comfort during lengthy procedures. Adjustable straps with hook and loop fastening are provided for increased patient stability. Twin-wheel locking casters provide smooth transport and secure table positioning.

#### SPECIFICATIONS:

Dimensions: 82" l x 29.5" w x 35" h (209 x 75 x 89 cm)

Tabletop: 80" l x 24" w (203 x 61 cm)

Radiolucent Area: 75" l x 24" w (191 x 61 cm)

Tabletop Material: Carbon fiber

Mattress: Seamless, 2" thick (5 cm)

Attenuation:

Tabletop: 1 mm Aluminum equivalence

Mattress: .5 mm Aluminum equivalence

Wheels:

Front – 5" (12.7 cm) locking twin-wheel casters

Rear – 3" (7.62 cm) locking twin-wheel casters

Patient Restraints: two body straps

Finish: Powder coat

Patient Capacity: 500 lb (227 kg); weight tested to four times the patient load rating.

Weight: 125 lb (57 kg)

Warranty: Two years



**056-025** Table, Surgical C-Arm, Fixed Height

Optional:

**056-865** Arm Board, Carbon Fiber, one arm

**056-866** Arm Board, Carbon Fiber, two arms

**058-864** Arm Board, Surgical, Carbon Fiber

*As required by law, all Biodex C-Arm tables are  
CDRH registered and FDA Device Listed.*